



Designation: D3895 – 19

# Standard Test Method for Oxidative-Induction Time of Polyolefins by Differential Scanning Calorimetry<sup>1</sup>

This standard is issued under the fixed designation D3895; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon ( $\epsilon$ ) indicates an editorial change since the last revision or reapproval.

## 1. Scope\*

1.1 This test method outlines a procedure for the determination of oxidative-induction time (OIT) of polymeric materials by differential scanning calorimetry (DSC). It is applicable to polyolefin resins that are in a fully stabilized/compounded form.

1.2 The values stated in SI units are to be regarded as the standard.

1.3 *This standard does not purport to address all of the safety concerns, if any, associated with its use. It is the responsibility of the user of this standard to establish appropriate safety and health practices and determine the applicability of regulatory limitations prior to use.* Specific hazards information is given in Section 8.

NOTE 1—This standard and ISO 11357–6 2013 address the same subject matter, but differ in technical content.

1.4 *This international standard was developed in accordance with internationally recognized principles on standardization established in the Decision on Principles for the Development of International Standards, Guides and Recommendations issued by the World Trade Organization Technical Barriers to Trade (TBT) Committee.*

## 2. Referenced Documents

2.1 *ASTM Standards:*<sup>2</sup>

D4703 Practice for Compression Molding Thermoplastic Materials into Test Specimens, Plaques, or Sheets

E473 Terminology Relating to Thermal Analysis and Rheology

E691 Practice for Conducting an Interlaboratory Study to Determine the Precision of a Test Method

E967 Test Method for Temperature Calibration of Differen-

tial Scanning Calorimeters and Differential Thermal Analyzers

E2935 Practice for Conducting Equivalence Testing in Laboratory Applications

## 3. Terminology

3.1 *Definitions*—Definitions of terms applying to thermal analysis appear in Terminology E473.

3.2 *Definitions of Terms Specific to This Standard:*

3.2.1 *oxidative induction time*—a relative measure of a material's resistance to oxidative decomposition; it is determined by the thermoanalytical measurement of the time interval to onset of exothermic oxidation of a material at a specified temperature in an oxygen atmosphere.

3.2.2 Abbreviations:

3.2.3 HDPE—high-density polyethylene.

3.2.4 LDPE—low-density polyethylene.

3.2.5 LLDPE—linear low-density polyethylene.

3.2.6 OIT—oxidative induction time.

## 4. Summary of Test Method

4.1 The sample to be tested and the reference material are heated at a constant rate in an inert gaseous environment (nitrogen). When the specified temperature has been reached, the atmosphere is changed to oxygen maintained at the same flow rate. The specimen is then held at constant temperature until the oxidative reaction is displayed on the thermal curve. The OIT is determined from the data recorded during the isothermal period. The time interval from when the oxygen flow is first initiated to the oxidative reaction is referred to as the induction period.

4.1.1 The end of the induction period is signaled by an abrupt increase in the specimen's evolved heat or temperature and will be recorded as an exothermic event by a differential scanning calorimeter (DSC).

4.2 The type of containment system used depends on the intended application use of the material being tested. Polyolefins used in the wire and cable industry typically require copper or aluminum pans, whereas polyolefins used in geomembrane and vapor-barrier film applications exclusively use aluminum pans.

<sup>1</sup> This test method is under the jurisdiction of ASTM Committee D20 on Plastics and is the direct responsibility of Subcommittee D20.30 on Thermal Properties.

Current edition approved May 1, 2019. Published June 2019. Originally approved in 1980. Last previous edition approved in 2014 as D3895 – 14. DOI: 10.1520/D3895-19.

<sup>2</sup> For referenced ASTM standards, visit the ASTM website, www.astm.org, or contact ASTM Customer Service at service@astm.org. For *Annual Book of ASTM Standards* volume information, refer to the standard's Document Summary page on the ASTM website.

\*A Summary of Changes section appears at the end of this standard

4.3 Unless otherwise specified, the analysis temperature used in this test is set at 200.0°C. For samples that have relatively low or high stabilization levels, it is permissible to adjust the temperature (typically between 180 and 220°C) to yield a thermal curve that can be interpreted and analyzed easily.

## 5. Significance and Use

5.1 The OIT is a qualitative assessment of the level (or degree) of stabilization of the material tested. This test has the potential to be used as a quality control measure to monitor the stabilization level in formulated resin as received from a supplier, prior to extrusion.

NOTE 2—The OIT measurement is an accelerated thermal-aging test and as such can be misleading. Caution should be exercised in data interpretation since oxidation reaction kinetics are a function of temperature and the inherent properties of the additives contained in the sample. For example, OIT results are often used to select optimum resin formulations. Volatile antioxidants may generate poor OIT results even though they may perform adequately at the intended use temperature of the finished product.

NOTE 3—There is no accepted sampling procedure, nor have any definitive relationships been established for comparing OIT values on field samples to those on unused products, hence the use of such values for determining life expectancy is uncertain and subjective.

## 6. Apparatus

6.1 *Differential Scanning Calorimeter*—As a minimum requirement, the thermal analysis equipment shall be capable of measuring heat flow of at least 10-mW full scale. The instrument recorder shall be capable of displaying heat flow or temperature differential on the *Y*-axis and time on the *X*-axis. The time base must be accurate to  $\pm 1\%$  and be readable to 0.1 min.

NOTE 4—The OIT test is a function of a particular compound's stabilizer system and should not be used as a basis of comparison between formulations that might contain, different resins, stabilizers, or additive packages, or all of these.

6.2 *Gas-Selector Switch and Regulators*, for high-purity nitrogen and oxygen. The distance between the gas-switching point and the instrument cell shall be such that the time required to transition to an oxygen environment is less than one minute. At a flow rate of 50 mL/min, this equates to a maximum switching volume of less than 50 mL.

6.3 *Analytical Balance*, 0.1-mg sensitivity.

6.4 *Bore-Hole Cutter*, 6.4-mm diameter.

6.5 *Specimen-Encapsulating Press*.

6.6 *Forceps, Scalpel, and Cutting Board*.

6.7 *Electronic Mass Flow Controller, Rotometer (Calibrated) or Soap-Film Flowmeter*, for gas-flow calibration.

6.8 *Specimen or Sample Pans*—Degreased aluminum or oxidized-copper pans (6.0 to 7.0-mm diameter, 1.5-mm height). Use the appropriate pan type for the material being tested.

NOTE 5—Aluminum lids are required for temperature calibration.

NOTE 6—Both types of pans are commercially available. Alternatively, the copper pans can be fabricated manually. Details on copper pan preparation and oxidation as well as instructions for aluminum pan

conditioning (degreasing) are given in [Annex A2 – Annex A4](#).

NOTE 7—The material composition of the specimen holder can influence the OIT test result significantly (that is, including any associated catalytic effects).

6.9 *Compression-Molding Device with Heated Platens*.

6.10 *Spacer Plates, Shim Stock, Caul Plates, etc.*

6.11 *Polyethylene Terephthalate Film (PET) or Polytetrafluoroethylene (PTFE) Coated Cloth*, for sample-plaque preparation.

6.12 *Thickness Gauge*.

6.13 *Laboratory Gas Burner*, for copper-pan oxidation.

6.14 *Boiling Flask*, with condenser and heating mantle.

6.15 *Forced-Air Oven*.

## 7. Reagents and Materials

7.1 All chemical reagents used in this procedure shall be analytical grade unless otherwise specified.

7.2 *Oxygen*—Ultra-high-purity grade (extra dry).

7.3 *Nitrogen*—Ultra-high-purity grade (extra dry).

7.4 *Aluminum Pan Degreasing Solvent*.

7.5 *Indium* (99.999 % purity).

7.6 *Tin* (99.999 % purity).

## 8. Hazards

8.1 Oxygen is a strong oxidizer that accelerates combustion vigorously. Keep oil and grease away from equipment using or containing oxygen.

8.2 The use of pressurized gas requires safe and proper handling.

## 9. Sampling

9.1 The following sample preparation procedures are recommended: the test sample is compression molded into sheet format (thickness of  $250 \pm 15 \mu\text{m}$ ) prior to analysis to yield consistent sample morphology and weight. Specimen disks cut from the sheet using 6.4-mm diameter bore-hole cutter will have a weight of approximately 5 to 10 mg, depending on sample density.

NOTE 8—If the sample requires homogenization prior to analysis, the melt compounding procedure given in [Appendix X1](#) is recommended. Poor sample uniformity will affect test precision adversely.

9.1.1 Obtain the required mass of the sample and place the material in the center of the appropriately sized spacer between two sheets of Polyethylene Terephthalate or PTFE coated cloth and two caul plates.

9.1.2 Place the assembly into the compression-molding device. The preheat and pressing temperature is 160°C for polyethylene and 190°C for polypropylene.

9.1.3 Heat the sample with appropriate pressure and time settings to obtain a plaque with uniform thickness.

9.1.4 Remove the plaque assembly and place it between two thick steel plates (heat sink) and cool the plaque to ambient temperature. Alternatively, quench the plaque in ice water.

9.1.5 Determine the average thickness of the sheet to ensure that it is within the allowable tolerances.

Sample: Tin calibration  
 Size: 5.1700 mg  
 Method: TIN CALIBRATION  
 Comment: 10°C/min to 220°C, then 1°C/min to 240°C

DSC

File: TT-TIN. 02  
 Operator: VD  
 Run Date: 2–Nov–90 13: 03

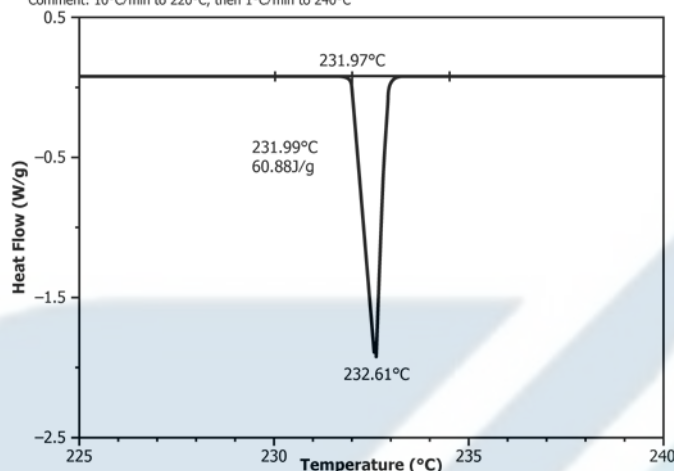


FIG. 1 Indium and Tin Melting Thermal Curves

9.1.6 Use the bore-hole cutter to punch out a disk from the plaque and record the specimen weight.

9.1.7 Place the specimen disk into the appropriate pan type. Use an identical empty pan as the reference. (Do not crimp or seal the pans.)

NOTE 9—If controlled cooling is not necessary, the option to prepare the test sample using Practice D4703, Annex 1 (Procedure C), is acceptable.

## 10. Procedure

10.1 *Instrumental Calibration*—This procedure uses a two-point calibration step. Indium and tin are used as the calibrants since their respective melting points encompass the specified analysis temperature range (180 to 220°C). Calibrate the instrument in accordance with the manufacturer’s instructions using the following procedure. Calibrate the instrument at least once per month or before use if longer than one month.

10.1.1 Place  $5 \pm 0.5$  mg of indium/tin into an aluminum sample pan. Place an aluminum cover over the pan, and seal using the encapsulating press. Prepare an empty sealed pan to be used as the reference. Place the specimen and reference pans into their respective locations in the instrument cell.

10.1.2 Turn on the nitrogen-gas flow at a rate of 50 mL/min (with an absolute pressure of 140 kPa).

10.1.3 Use the following melting profiles:

*Indium:* = ambient to 145°C at 10°C/min, 145 to 165°C at 1°C/min

*Tin:* = ambient to 220°C at 10°C/min, 220 to 240°C at 1°C/min

NOTE 10—The specified heating rates are for calibration use only.

10.1.4 Adjust the temperature-calibration software (or potentiometer) to set the melting point at 156.63 and 231.97°C (see Practice E967) for indium and tin, respectively. The melting point of the calibrant is defined as the intercept of the

extended baseline and the extended tangent to the first slope of the endotherm (that is, the onset). See Fig. 1.

NOTE 11—An inadequate melting thermal curve is occasionally obtained due to poor surface contact of the calibrant material to the pan surface. If this occurs, repeat the calibration step. (After one melting/crystallization cycle the calibrant material should coat the bottom of the pan evenly.)

### 10.2 Instrument Operation:

10.2.1 Load the specimen and reference pans into the cell.

10.2.2 Allow 5 min for a nitrogen prepurge prior to beginning the heating cycle to eliminate any residual oxygen. Commence programmed heating of the specimen (under nitrogen flow of  $50 \pm 5$  mL/min) from ambient temperature to 200°C (set point) at a rate of 20°C/min.

10.2.3 When the set temperature has been reached, discontinue programmed heating and equilibrate the sample for 5 min at the set temperature. Turn on the recorder. If the instrument being used does not have an isothermal temperature-control-mode feature, follow the alternate procedure outlined in Annex A1 or alternatively ensure accurate temperature control by monitoring and adjusting continually, as required.

10.2.4 Once the equilibrium time has expired, change the gas to oxygen at a flow rate of  $50 \pm$  mL/min. (Record this event.) This changeover point to oxygen flow is considered the zero time of the experiment.

10.2.5 Continue isothermal operation until the maximum exotherm has been reached to allow a complete examination of the entire exotherm. (see Fig. 2 and Fig. 3). At the tester’s discretion, it is acceptable to terminate the test at a predetermined heat flow change provided that data are available to support the alternative. It is also acceptable to terminate the test if time requirements stated in the product’s specification have been met.

10.2.6 Upon completion of the test, switch the gas selector back to nitrogen and cool the instrument to ambient temperature. If additional testing is being conducted, cooling the



Sample: Indium calibration  
 Size: 5.1700 mg  
 Method: INDIUM CALIBRATION  
 Comment: 10°C/min to 145°C, then 1°C/min to 165°C

**DSC**

File: TTINDIUM.02  
 Operator: VD  
 Run Date: 1-Nov-90 15:44

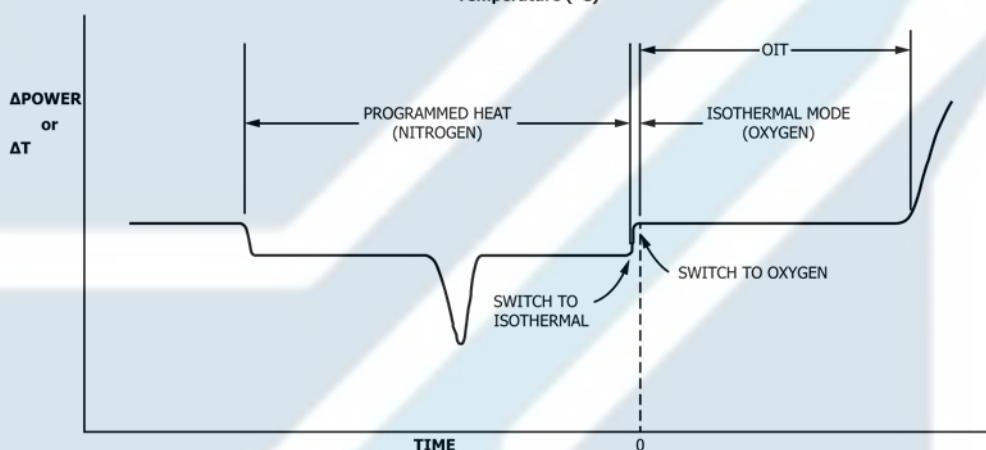
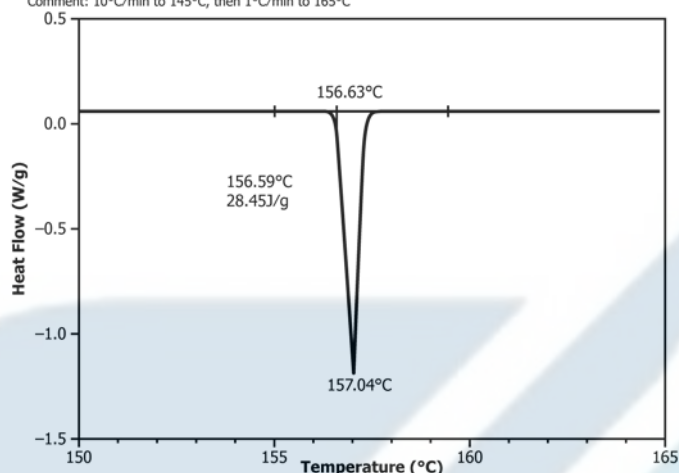


FIG. 2 Evaluation of OIT from Recorded-Time-Based Thermal Curve

Sample: OIT-2  
 Size: 7.0000 mg  
 Method: OIT at 200°C  
 Comment: 02, 50 ML/MIN, 7 PSI

**DSC**

File: A: OIT-2.05  
 Operator: SB  
 Run Date: 7-Jun-91 11:37

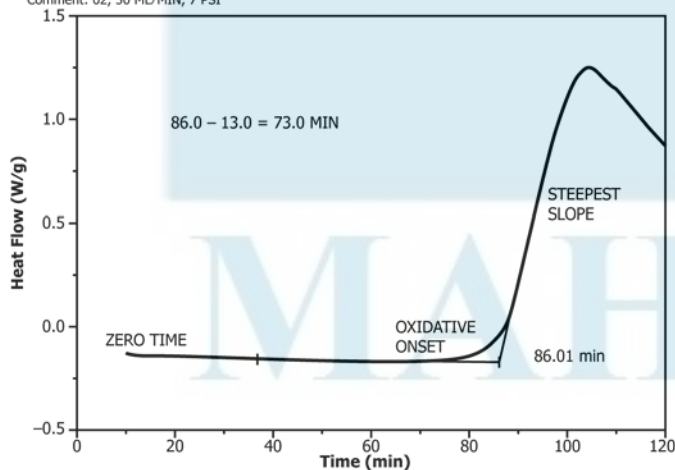


FIG. 3 Determination of OIT

10.2.7 Test frequency is established by the user. As a minimum requirement, samples are tested in duplicate with the mean value reported.

10.2.8 Clean the DSC cell of contamination by heating to 500°C for 5 min in air (or oxygen) prior to conducting measurements and between the testing of different formulations.

10.3 *Thermal Curve Analysis*—The data is plotted with the heat-flow signal normalized to sample mass (that is, W/g) on the y-axis, versus time on the x-axis. Expand the x-axis as much as possible to facilitate analysis.

10.3.1 Extend the recorded baseline beyond the oxidative reactive exotherm. Extrapolate the steepest linear slope of this exotherm to intercept the extended baseline (see Fig. 3).

10.3.2 The OIT is measured to within ±0.1 min from zero time to the intercept point.

10.3.3 The tangent method used to measure the oxidation time is the preferred method, but the selection of the appropriate tangent to the exotherm sloped line is, at times, difficult if the exothermic peak has a leading edge. It is possible that exothermic peaks with leading edges occur if the oxidation reaction is slow.

NOTE 12—The oxidation thermogram is indicative of the additives used

instrument cell below 60 to 70°C is sufficient to avoid any premature thermal oxidation of the sample.

**TABLE 1 Precision and Repeatability—Oxidative Induction Time (min) of Polyethylenes**

Material	Average	$S_r$	$S_R$	$r$	$R$
LDPE	24.0	2.8	4.3	8.0	12.2
LDPE	83.4	9.2	17.4	25.7	48.8
LLDPE	120	7.8	14.6	21.9	40.8
HDPE	163	8.3	21.7	23.2	60.9

in the polymer. Thus, oxidation may not occur as a smooth transition and multiple steps or slopes may be produced. When this occurs, resampling and retesting is recommended to ensure that the testing is representative of the oxidation process.

NOTE 13—If multiple slopes result from the oxidation process, OIT needs to be defined to accurately reflect the oxidation of the polymer. It is up to the user to determine which slope best represents the material property for an application. It must be noted in the report if the tangent line is not drawn from the first steep slope.

NOTE 14—If the selection of the appropriate baseline is not obvious using the tangent method, try the offset method. Draw a second baseline parallel to the first baseline at a distance of 0.05 W/g above the first baseline. The intersection of this second line with the exotherm signal is defined as the onset of oxidation. The time from this intersection to zero time is the OIT.

NOTE 15—Other gas-flow rates or specimen weights, or both, may be used if shown to be equivalent.

## 11. Report

11.1 Report the following information:

11.1.1 Identification of the sample,

11.1.2 Type of sample pan,

11.1.3 Test temperature used,

11.1.4 Measurement technique used (that is, tangent method or offset method, OIT not determined from the first steep slope), and,

11.1.5 Average OIT in minutes (three significant digits).

## 12. Precision and Bias<sup>3</sup>

12.1 The data in Table 1 is based on a round robin conducted in 1991 in accordance with Practice E691, involving four polyethylene materials tested by eleven laboratories. The test materials were obtained from three laboratories. The OIT of each sample was determined in duplicate on two separate days (that is, a total of four OIT determinations). Of the fourteen laboratories that participated in the round-robin exercise, eleven laboratories used heat flux-type instruments and three laboratories used power compensation-type instruments.

<sup>3</sup> Supporting data are available from ASTM Headquarters. Request RR:D20-1178.

NOTE 16—The precision and bias data for other polyolefins may differ from the polyethylene material data shown in Table 1.

NOTE 17—The sample preparation steps for each test material (that is, sheet compression molding) were performed by one laboratory. Each participating laboratory received a set of specimen disks die cut from each of the molded sheets and samples of indium and tin required for instrument calibration. Each participating laboratory furnished its own specimen holders.

12.2 The following explanations of  $r$  and  $R$  (12.3 – 12.3.3) are intended to present only a meaningful way of considering the approximate precision of this test method. Do not apply the data in Table 1 rigorously to acceptance or rejection of material, as those data are specific to the round robin, and not representative of other lots, condition, materials, or laboratories. Users of this test method are encouraged to apply the principles outlined in Practice E691 to generate data specific to their laboratory and materials, or between specific laboratories. The principles of 12.3 – 12.3.3 would then be valid for such data.

12.3 *Concept of  $r$  and  $R$* —If  $S_r$  and  $S_R$  have been calculated from a large enough body of data, and for test results that were averages from testing four specimens, the following applies:

12.3.1 *Repeatability*—The value below which the absolute difference between two individual test results obtained under repeatability conditions is likely to be expected to occur with a probability of approximately 0.95 (95 %). “ $r$ ” is the interval representing the critical difference between two test results for the same material, obtained by the same operator using the same equipment on the same day in the same laboratory.

12.3.2 *Reproducibility*—The value below which the absolute difference between two test results obtained under reproducibility conditions is likely to be expected to occur with a probability of approximately 0.95 (95 %). “ $R$ ” is the interval representing the critical difference between two test results for the same material, obtained by different operators using different equipment in different laboratories.

12.3.3 Conducting equivalence testing on numerical data from two sources shall be conducted in accordance with Practice E2935 or any known method for judging the equivalence of two means.

NOTE 18—Example, a t-test.

12.4 *Bias*—There are no recognized standards by which to estimate the bias of this test method.

## 13. Keywords

13.1 differential scanning calorimetry; oxidative induction time; polyolefins

(Mandatory Information)

A1. ALTERNATE OIT DETERMINATION

A1.1 If the DSC instrument being used does not have an isothermal temperature-control-mode feature, the following procedure is used for OIT determination. Commence programmed heating (under a nitrogen flow of  $50 \pm 5$  mL/min) from ambient temperature and record the thermal curve. Continue heating until  $200.0 \pm 2^\circ\text{C}$  is reached. Discontinue programmed heating and equilibrate the specimen for 5 min to a constant temperature. Measure and record this temperature to  $\pm 0.1^\circ\text{C}$ . Once the equilibrium time has expired, change the gas to oxygen at a flow rate of  $50 \pm 5$  mL/min. This changeover point to oxygen flow is considered the zero time of the experiment. Continue isothermal operation until at least 2 min have elapsed after the steepest point of the exotherm has been displayed (see Fig. 2). Upon completion of the test, switch the gas selector back to nitrogen and cool the instrument to ambient temperature. Repeat the entire test on a new sample, thereby generating a second exotherm. The isothermal portion of these two tests must be different from each other by at least  $1^\circ\text{C}$ , but within  $\pm 2^\circ\text{C}$  of the set temperature (for example,  $200.0 \pm 2^\circ\text{C}$ ).

A1.2 Extend the recorded baseline beyond the oxidative reactive exotherm. Extrapolate the slope of this exotherm to intercept the extended baseline (see Fig. 3). The OIT is measured to within  $\pm 0.1$  min from zero time to the intercept

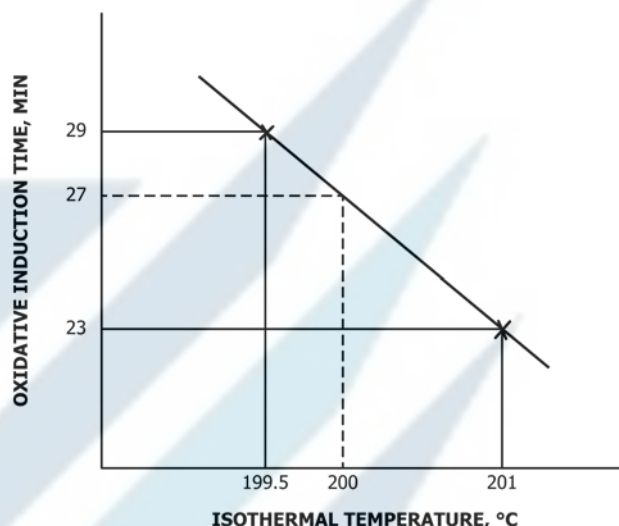


FIG. A1.1 Example of the Linear-Interpolation Technique Used to Estimate the OIT at  $200.0^\circ\text{C}$

point. Make a plot of OIT versus temperature. Interpolate between the two actual temperature values until an OIT is obtained that corresponds exactly to the set temperature (for example,  $200.0^\circ\text{C}$ ). A sample calculation is given in Fig. A1.1.

A2. PREPARATION OF COPPER PANS

A2.1 Punch out a 10-mm diameter disk from a 0.10-mm thick soft-copper strip. Using a punch and die tool, draw the copper disks to form shallow holders. The lower 1.5 to 1.8 mm of the side of the holder must be straight and even, but the edge

sections above will likely be irregular or ragged. The diameter of the holder at the base must be from 6.6 to 6.9 mm. Trim off the ragged edges.

A3. OXIDATION OF COPPER PANS

A3.1 Prepare fresh copper oxide coatings on the copper pans by placing each pan in the flame of a gas burner until it begins to glow (approximately 5 to 10 s). Rotate the pan in the flame slowly to ensure a uniform oxide coating. Remove the pan from the flame and allow it to cool to room temperature (approximately 30 s).

#### A4. ALUMINUM PAN DEGREASING

A4.1 Condition aluminum pans and lids (lids used for calibration only) prior to use. Reflux the pans and lids in a low-boiling solvent such as cyclohexane for 30 min, and then dry them in a forced-air oven at 40 to 50°C.

### APPENDIX

#### (Nonmandatory Information)

#### X1. SAMPLE HOMOGENIZATION

X1.1 If desired, the test sample may be homogenized prior to OIT analysis to minimize local variations of stabilizer concentration. A torque rheometer equipped with a mixing head is suitable for such purposes.

X1.2 Each sample should be blended at 60 r/min for 10 min, using a chamber temperature of 150 to 160°C for polyethylene resins and 180 to 190°C for polypropylene resins. The mixing

chamber should be purged with nitrogen to suppress polymer and additive degradation.

X1.3 After a mixing time of 10 min, quench-cool the samples in liquid nitrogen or ice water. After sufficient cooling, compression mold the required amount of material into sheet format.

### SUMMARY OF CHANGES

Committee D20 has identified the location of selected changes to this standard since the last issue (D3895 - 14) that may impact the use of this standard. (May 1, 2019)

- |   |   |
|---|---|
| <p>(1) Revised 4.3 for clearer understanding of the required test temperature and 9.1.</p> <p>(2) Modified 12.3.1, 12.3.2, and 12.3.3 for definitions of repeatability and reproducibility.</p> | <p>(3) Added <b>Note 18</b>.</p> <p>(4) Added reference Practice <b>E2935</b> to <b>2.1</b>.</p> <p>(5) Adjusted reference footnotes.</p> |
|---|---|

*ASTM International takes no position respecting the validity of any patent rights asserted in connection with any item mentioned in this standard. Users of this standard are expressly advised that determination of the validity of any such patent rights, and the risk of infringement of such rights, are entirely their own responsibility.*

*This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and if not revised, either reapproved or withdrawn. Your comments are invited either for revision of this standard or for additional standards and should be addressed to ASTM International Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical committee, which you may attend. If you feel that your comments have not received a fair hearing you should make your views known to the ASTM Committee on Standards, at the address shown below.*

*This standard is copyrighted by ASTM International, 100 Barr Harbor Drive, PO Box C700, West Conshohocken, PA 19428-2959, United States. Individual reprints (single or multiple copies) of this standard may be obtained by contacting ASTM at the above address or at 610-832-9585 (phone), 610-832-9555 (fax), or service@astm.org (e-mail); or through the ASTM website (www.astm.org). Permission rights to photocopy the standard may also be secured from the Copyright Clearance Center, 222 Rosewood Drive, Danvers, MA 01923, Tel: (978) 646-2600; http://www.copyright.com/*